

Isolator Glove Sleeve Clamps


CRL Part # HDSC



At left – a HDPE Isolator Glove Sleeve Clamp used to replace gloves while an isolator is in operation. The device is applied to the glove sleeve squeezing the lumen of the sleeve thereby separating the glove and a portion of the sleeve from the rest of the isolator providing a microbiological barrier during the glove changing process.

At right – Dimensions:

Height: 2.50 inches
Length: 23 inches
Width: 1 inch
Weight: 3 lbs

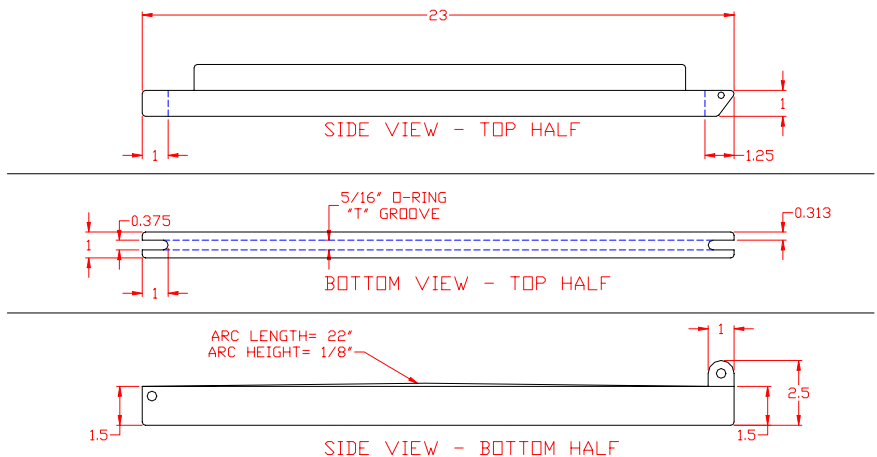

Isolator Glove Sleeve Clamp

CRL Part # HDSC